

NEW

astris-pme.com.au 1300 131 884

13

Stander Prone/Supine/Multi-Position

It's finally here – the Stander that has it all. We're thrilled to introduce a total redesign of our Supine and Prone Stander families. This new Stander can be configured as supine, prone or multi-position.

All models allow hip abduction up to 30° to promote hip joint development. Independent adjustments for hip and knee flexion and leg length accommodate contractures and discrepancies. A broad range of positioning components allow users to remain comfortably upright and active for longer. The frame can be raised and lowered to position a child at eye level with peers, or in easy reach of activities.

All this function is contained within an elegant and approachable design that centers the child while opening worlds of possibility.

Find important details and studies on the therapeutic benefits of standing at: www.rifton.com/standers

ability to fully customize joint positioning through the lower extremities for maximal comfort and alignment is a game changer! I have a student who only tolerated five to ten minutes of static standing in a less adjustable stander. When we tried the new Rifton Stander the student was happy and comfortable for nearly an hour, and we only got him out because we ran out of time.

Linda K. Miller, PT, DPT MOVE International Trainer Prospect Center, NY

•••••• Size 2 Prone Stander configuration (S420)

Stander features

Any of the new Standers can be easily switched between prone and supine.



Grifton

Multi-Position Stander



The Stander can be purchased in the prone, supine or multi-position configuration. If the multi-position configuration is ordered, the conversion between prone and supine is quick, easy and tool free.

Quick and easy conversion between prone and supine



- **1** Add/Remove the upper trunk support.
- **2** Switch the tray direction.
- Swap the supine hip strap and prone pelvic harness.
- **4 Add/Remove** the supine knee straps.
- **5** Turn the sandals 180°.



Quick, precise positioning



Advanced legs have independent leg length, hip angle and knee angle adjustments.



Up to 30° of abduction.



Tilt adjustments from 0°-90°.



Hip abduction

Independent leg abduction, up to 30° in each leg, is now standard. It's what clinicians everywhere asked us for to promote healthy hip development.

G rifton.



The size 2 and 3 Stander's advanced legs have independent hip and knee angle adjustment to accommodate contractures. Separate upper and lower leg length adjustments ensure that users with leg length discrepancies are comfortably supported.





Move between eye-level and activity level within seconds.

- ▲ When the frame is lowered, users can interact at eye-level with their peers.
- With the frame raised, users can access activities and environments.

Grifton.

Quick, safe, easy transfers

Prone

Sit-to-stand transfers are easy in the prone configuration.



Supine





Pads and supports swing away for easy lateral transfers in the supine configuration.

Components



Knee cuffs (Prone only)

The knee cuffs provide lateral stabilization and rotational control and prevent knee hyperextension.







The upper trunk support provides security and safety while still allowing freedom of arm movement. It is uniquely designed to comfortably support all body types, including female users.

Note: Available for size 3 only.

Head laterals (Supine only) The head laterals mount on the headrest to provide lateral head support. They adjust independently in height and width.










The butterfly harness can be used instead of the trunk strap for greater upper trunk stabilization. The slim cut butterfly harness (size 3 only) is an alternative especially designed for larger female clients.





Rotation control

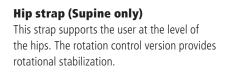





Soft strap



Knee supports (Supine only)

The soft knee straps form around the knee and have a cutout for the patella. For more rigid support, choose the adjustable split knee supports which swing away for easy loading.





STANDER



Trunk strap

The trunk strap supports the user's chest. The rotation control version wraps around the user and provides rotational stabilization.

Leg types

Standard legs allow for hip abduction, and independent leg length adjustments. Advanced legs are optional on sizes 2 and 3 and include hip and knee flexion in both prone and supine configurations.

Advanced





Tray

The tray can be used in prone and supine configuration. Its position and angle adjust easily, and it swings away for transfer.

Additional trunk laterals

The laterals attach at chest or hip level and swing away for transfer. The pads rotate for individual height adjustment.





Medial/lateral thigh supports The medial/lateral thigh supports provide abduction or adduction. The mounting bracket is off-center, so the pad has two possible height positions depending which side it is installed on.





Sandal lift kit The smallest users ma

The smallest users may require the S409 sandal lift kit for best positioning. **Note:** Available for size 1 only.









Stander dimensions

User dimensions (inches)	S410 Size 1	S420 Size 2	S430 Size 3						
Height	28-42*	38-53	48-65						
Key user dimension: height Select the appropriate stander by the user's overall height. Choose the model that allows for growth. *The smallest users may require the S409 Sandal lift kit for best positioning. Important: User's weight must not exceed the maximum working load.									
Item dimensions (inches)	S410 Size 1	S420 Size 2	S430 Size 3						
Base length x width	32½ x 21½	39 x 27	43 x 29						
Height of horizontal body support above floor	22-28	221/2-361/2	231/2-38						
Sandal to knee height	6¾-11	10-151/4	121/4-19						
Knee to hip height	61/4-101/2	81/2-131/2	12-171/4						
Sandal to top of headrest	311/2-43	41–55	49-661/2						
Width between trunk laterals	51/2-10	7-111/2	9-14						
Width between medial/lateral thigh supports	21/2-4	31⁄4-5	4-51/2						
Basic item weight PRONE (lb)	47	72	104						
Basic item weight SUPINE (lb)	52	78	115						
Max. working load (lb)	75	100	185						



New! Size 3 Stander at its smallest and largest settings

Grifton.

STANDER

2025 order form

Use dimension chart to select appropriate size.



Scan for sizing and pricing:

ose dimension chart to select appropriate size.			S410 size 1		S420 size 2			S430 size 3				
* Fi	ame				S411	\$1,940 🗌	S421	\$2,550		S431	\$3,160	
			£	Blue		Blue			Blue			
				CHOOSE	Red		Red	••••		Red	••••	
	Prone (with pelvic harness)				S412	\$305	S422	\$315		S432	\$325	
		New! Prone upper	trui	nk support (recommended)	N/A		N/A			S408	\$250	
		*Trunk strap ♀	£	Standard	S462	\$305 🗌	S463	\$315		S464	\$325	
			CHOOSE	Rotation control	S466	\$305	S467	\$315		S468	\$325	
		Prone knee cuffs (pa			S485	\$205	S485	\$205		S486	\$225	
	Supine (with upper				S413	\$355	S423	\$365		S433	\$375	
	trunk support)	Head laterals (pair)	₽	Small	S497	\$205	S497	\$205		N/A		_
			CHOOSE	Large	N/A	, L	\$498	\$215		S498	\$215	 _
		* Trunk straps		Butterfly harness (standard)	\$456	\$305 🗌	\$457	\$315		\$458	\$325	
		CHOOSE	ç	Butterfly harness (slim-cut)	N/A	,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,	N/A	4515		S459	\$325	
			sool	Trunk strap (standard)	\$462	\$305 🗌	S463	\$315		S464	\$325	
*			m	Trunk strap (rotation control)	S466	\$305	S467	\$315		S468	\$325	
CHOOSE		* Hip strap	c	Standard	S472	\$305	S473	\$315		S474	\$325	
			CHOOSE	Rotation control	S472	\$305	S477	\$315		S478	\$325	
CONFIGURATION			-		S401	\$205	S402			S403		
FIG		kilce supports	CHOOSE	Soft straps (pair)		\$205		\$215 ¢42E			\$225	
JRA			m	Split supports (pair)	N/A		S487	\$435		S488	\$460	
	Multi-Position (with head support, pelvic harness)	Hood laterals	0	C II	S414	\$665	S424	\$685		S434	\$705	
-		Butterfly harness * Trunk strap * Hip strap * Knee supports	CHOOSE	Small	S497	\$205 🗌	S497	\$205		N/A	****	
				Large	N/A		S498	\$215		S498	\$215	
			CHOOSE	Standard	S456	\$305 🗌	S457	\$315		S458	\$325	
				Slim-cut	N/A		N/A			S459	\$325	
			CHOOSE	Standard	S462	\$305 🗌	S463	\$315		S464	\$325	
				Rotation control	S466	\$305 🗌	S467	\$315		S468	\$325	
			CHOOSE	Standard	\$472	\$305 🗌	S473	\$315		S474	\$325	
				Rotation control	S476	\$305 🗌	S477	\$315		S478	\$325	
			CHOOSE	Soft straps (pair)	S401	\$205 🗌	S402	\$215		S403	\$225	
				Split supports (pair)	N/A		S487	\$435		S488	\$460	
				nk support (recommended)	N/A		N/A			S408	\$250	
	Prone knee cuffs (pair)		S485	\$205 🗌	S485	\$205		S486	\$225			
* Trunk laterals Small		Small	S494	\$205 🗌	S494	\$205		N/A	.			
			OSE	Large	N/A		S495	\$215		S495	\$215	
* Leg type		CHOOSE	Standard	S415	\$665 🗌	S425	\$685		S435	\$705	C	
			OSE	Advanced	N/A		S426	\$990		S436	\$1,010	Ľ
Tray		S491	\$560 🗌	S491	\$560		S492	\$610				
Additional trunk laterals		сно	Small	S494	\$205 🗌	S494	\$205		N/A			
			CHOOSE	Large	N/A		S495	\$215		S495	\$215	Ľ
Medial/lateral thigh supports Small (pair) Large		Small	S482	\$205 🗌	S482	\$205		N/A				
		OSE	Large	N/A		S483	\$215		S483	\$215		
Additional medial/lateral thigh		Small	S482	\$205 🗌	S482	\$205		N/A				
		DOSE	Large	N/A		S483	\$215		S483	\$215		
Size	Size 1 sandal lift kit			S409	\$55 🗌	NA			NA			
Headrest adapter plate kit			N/A		\$427	\$255		S427	\$255			
					1		5.21	4200		5.21	4200	

* Required

